

# 1650™

Manual operating table

---

Special Edition

For general surgery



OPERATING ROOM

DAY SURGERY

## Manual operating table

Robust and economic, the 1650 operating table has been one of Merivaara's best-selling products over the years. Now it has been made even more durable with a more stable base and new column design. It accommodates also heavier patients, and offers better access to the operated site thanks to the improved adjustment options.

Operating table 1650 is an economic and durable option for general surgery needs.

### Optional accessories

- Divided leg sections
- Narrow/special head rest
- Any accessories that can be fitted on the side rail, such as stirrups, special head rest, arm rests, etc.

### Adjustments

- Flex
- Reflex
- Trendelenburg
- Reverse Trendelenburg
- Hydraulic height adjustment
- Gas-spring assisted lateral tilt



General surgery, ENT/Struma



General surgery in lateral position, kidney/ thorax surgery with flex or separate kidney elevator



Laparoscopic surgery



Preparation for lithotomy position

### Lots of benefits for the surgical team

- Unobstructed working access
- Easy set up - faster turnaround times
- C-arm and X-ray compatible

Compared with the previous version of 1650, the new 1650 offers the following benefits:

- Higher safe working load (SWL)
- Enhanced lateral tilt for better patient positioning
- Better stability with new base and column design

For detailed specifications, see instruction manual and [www.merivaara.com](http://www.merivaara.com). You'll also find information on accessories and possible updates on the website.

### Technical specifications

Table top	4-5
Table top width (without rails)	500 mm
Weight	135 kg
Safe working load (SWL)	180 kg
Maximum lifting capacity	400 kg

#### ADJUSTMENTS RANGES

Height	730 mm – 1040 mm
Lateral tilt	± 20°
Trendelenburg	25°
Reverse Trendelenburg	18°
Back section	-40°...+70°
Leg section	-90°...+30°

#### ADJUSTMENTS

Back and leg section, lateral tilt and Trendelenburg	Mechanical, gas-spring assisted
Height	Hydraulic

- EN 60601-1, EN 60601-2-46 (when applicable)
  - 1650 operating table is classified into the class I according to European Council; Directive 93/42/EEC Class I (Annex IX)
- CE


# Solutions for efficient physical patient flow and well-being

Merivaara provides a wide range of hospital-grade furniture such as operating tables, medical lights, examination tables, trolleys and stretchers for transportation and day surgery, as well as birthing and patient beds. Additionally, our new concepts include both revolutionary integrated management systems and after-sales services.

Our versatile solutions and product offering enable convenient and cost-effective physical patient flow within the hospital, in day surgery and health clinics, as well as in nursing homes and home care. Our priority is to enhance comfort and well-being of both patients and their caregivers.

Established in 1901, Merivaara has more than 100 years of experience in designing and manufacturing hospital-grade furniture. Today our products are highly appreciated by users in more than 120 countries thanks to their ease-of-use, durability and ergonomic design.

Merivaara's medical products are CE-labelled and the company's quality management system complies with EU directives for medical device. It is certified by DNV according to ISO 9001:2008 and ISO 13485:2003 standards. The company's environmental management system is ISO 14001:2004 certified.

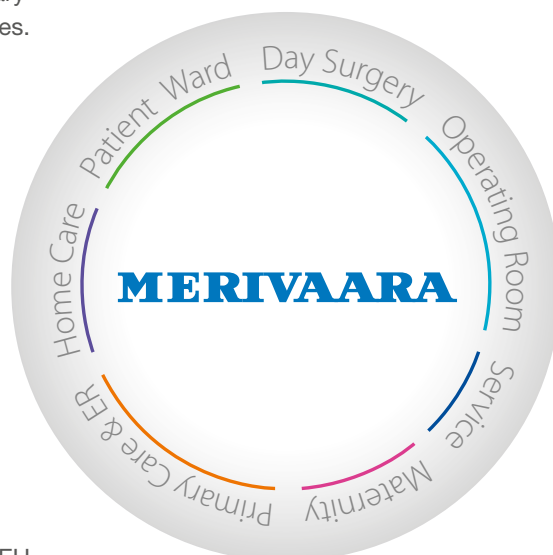

**Merivaara Corp.**  
Puustellintie 2, FI-15150 Lahti, FINLAND  
Tel. +358 3 3394 611, fax +358 3 3394 6144  
merivaara@merivaara.com  
www.merivaara.com

The status of the information, specifications and illustrations in this publication is indicated by the date of publication. Whilst the greatest care has been taken over its compilation, Merivaara Corp. is not responsible for any errors or omissions. Merivaara Corp. reserves the right to make changes to the technology, features, specifications and design of the equipment and models without notice. Please check the Merivaara website for the latest updated specifications and other information.

All trademarks are Merivaara property unless otherwise stated. Patents granted or pending.